

K-ASSAY®



KAMIYA BIOMEDICAL COMPANY

Rat Microalbumin ELISA

**For the quantitative determination of microalbumin
in rat biological samples**

Cat. No. KT-354

For research use only.

PRODUCT INFORMATION

Rat Microalbumin ELISA

Cat. No. KT-354

INTENDED USE

The Rat Microalbumin ELISA is a highly sensitive two-site enzyme-linked immunoassay (ELISA) for the quantitative determination of microalbumin in rat biological samples. **For research use only.**


INTRODUCTION

Albumin is an amazing polyfunctional protein contributing to homeostasis through mechanisms of hemodynamics, transport and nutrition. Albumin is found both intra- and extra-vascularly in all mammals and many lower vertebrates. It is a molecule of about 67,000 Daltons, synthesized by the liver. Normally only very trace amounts of albumin escape reabsorption by kidney glomeruli and are excreted into the urine which is then referred to as microalbumin. Many occult diseases can cause kidney damage which may result in excessive amounts of serum proteins, including microalbumin, to be excreted by the kidney and into the urine. This ELISA can be used to measure microalbumin in serum, tissue extracts and other biological fluids.

PRINCIPLE

The principle of the double antibody sandwich ELISA is represented in Figure 1. In this assay the microalbumin present in the sample reacts with the anti-microalbumin antibody, which has been absorbed to the surface of polystyrene microtiter wells. After the removal of unbound proteins by washing, anti-microalbumin antibody conjugated with horseradish peroxidase (HRP), are added. These HRP-conjugated antibody forms a complex with the previously bound microalbumin. Following another washing step, the enzyme bound to the immunosorbent is assayed by the addition of a chromogenic substrate, 3,3',5,5'-tetramethylbenzidine (TMB). The quantity of bound enzyme varies directly with the concentration of microalbumin in the sample tested; thus, the absorbance, at 450 nm, is a measure of the concentration of microalbumin in the test sample. The quantity of microalbumin in the test sample can be interpolated from the calibration curve constructed from the calibrators and corrected for sample dilution.

Figure 1. Principle of Assay

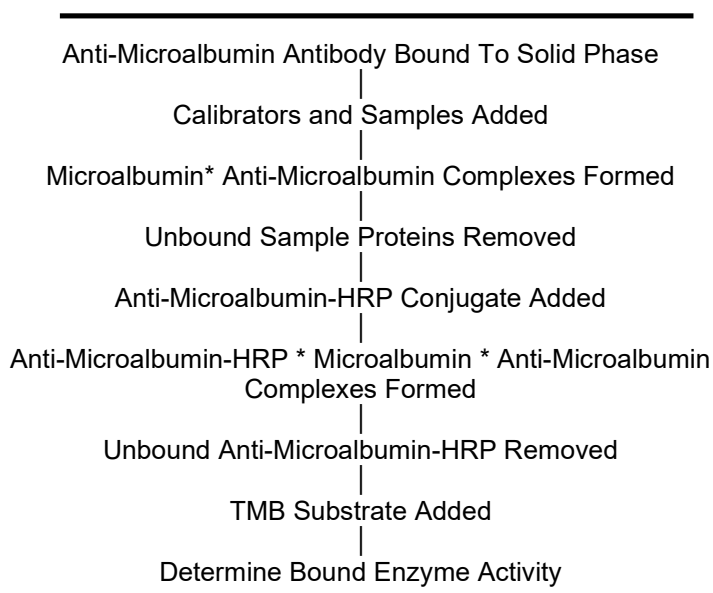

LIMITATION OF THE PROCEDURE

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC PURPOSES. IN VITRO USE ONLY.

Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the information contained in the package insert instructions and with adherence to good laboratory practice. Factors that might affect the performance of the assay include proper instrument function, cleanliness of glassware, quality of distilled or deionized water, and accuracy of reagent and sample pipettings, washing technique, incubation time or temperature. Do not mix or substitute reagents with those from other lots or sources.

COMPONENTS

1. **Diluent Concentrate**
One bottle containing 50 mL of a 5X concentrated diluent buffer.
2. **Wash Solution Concentrate**
One bottle containing 50 mL of a 20X concentrated wash solution.
3. **Enzyme-Antibody Conjugate Concentrate**
One vial containing 150 μ L of a 100X concentrated affinity-purified anti-rat microalbumin antibody conjugated with HRP in stabilizing buffer.
4. **TMB Substrate Solution**
One bottle containing 12 mL of TMB and hydrogen peroxide in citric acid buffer at pH 3.3.
5. **Stop Solution**
One bottle containing 12 mL of 0.3 M sulfuric acid.
WARNING: Avoid contact with skin.
6. **Microtiter Plate**
Twelve removable eight-well strips in a well holder frame. Wells are coated with affinity-purified anti-rat microalbumin antibody
7. **Rat Microalbumin Calibrator**
One vial containing a lyophilized Rat Microalbumin Calibrator.
8. **Positive Control**
One vial containing 50 μ L of serum with \leq 0.1% sodium azide. See the **Control Certificate** for the concentration.

MATERIALS REQUIRED BUT NOT PROVIDED

- Precision pipettes (2 μ L to 100 μ L) for making and dispensing dilutions
- Test tubes
- Squirt bottle or microplate washer/aspirator
- Distilled or de-ionized H₂O
- Microplate reader
- Assorted glassware for the preparation of reagents and buffer solutions
- Timer
- Centrifuge for sample collection
- Anticoagulant for plasma collection

STORAGE AND STABILITY

1. **Complete Kit**
The expiration date for the kit is stated on the outer label. The recommended storage temperature is 4°C. All components should be stable up to the expiration date if stored and used per this kit protocol insert. **Note: See long-term storage recommendations below for the Rat Microalbumin Calibrator and Positive Control**

1. **Diluent Concentrate**
The 5X Diluent Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions should be stored at 4°C.
2. **Wash Solution Concentrate**
The 20X Wash Solution Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions should be stored at 4°C.
3. **Enzyme-Antibody Conjugate**
Undiluted anti-Microalbumin-HRP conjugate should be stored at 4°C in the dark (protect from light) and **diluted immediately prior to use**. The undiluted conjugate solution is stable until the expiration date.
2. **TMB Substrate Solution**
The TMB Substrate Solution should be stored at 4°C in the dark (protect from light) and is stable until the expiration date.
4. **Stop Solution**
The Stop Solution should be stored at 4°C and is stable until the expiration date.
5. **Microtiter Plate**
Anti-rat microalbumin coated wells are stable until the expiration date of the kit and should be stored at 4°C in the sealed foil pouch with a desiccant pack.
6. **Rat Microalbumin Calibrator**
The lyophilized Rat Microalbumin Calibrator should be stored at 4°C until reconstituted. The reconstituted calibrator should be aliquoted and stored frozen. Avoid multiple freeze/thaw cycles. The working calibrator solutions should be prepared immediately prior to use.
7. **Positive Control**
For storage longer than 7 days, keep frozen until the expiration date. Storage for less than 7 days can be at 4°C. Avoid multiple freeze/thaw cycles. It is not recommended to store this product in a frost-free freezer.

PRECAUTIONS

1. Read the instructions carefully before beginning the assay.
2. This kit is for research use only.
3. Great care has been taken to ensure the quality and reliability of this product. However, it is possible that in certain cases, unusual results may be obtained due to high levels of interfering factors.
4. Preservatives
 - Positive Control contains $\leq 0.1\%$ sodium azide.
5. Azide and thimerosal at concentrations higher than 0.1% inhibit the enzyme reaction.
6. Other precautions:
 - Do not interchange kit components from different lots.
 - Do not use kit components beyond the expiration date.
 - Protect reagents from direct sunlight.
 - Do not pipette by mouth.
 - Do not eat, drink, smoke or apply cosmetics where reagents are used.
 - Avoid all contact with the reagents by using gloves.
 - Stop solution contains diluted sulfuric acid. Irritation to the eyes and skin is possible. Flush with water after contact.

SPECIMEN COLLECTION AND HANDLING

All blood components and biological materials should be handled as potentially hazardous. Follow universal precautions when handling and disposing.

If blood samples are clotted, grossly hemolyzed, lipemic, or the integrity of the sample is of concern, make a note and interpret results with caution.

The sample collection and storage condition listed below are intended as general guidelines. Sample stability has not been evaluated.

- **Serum samples** - Blood should be collected by venipuncture. The serum should be separated from the cells after clot formation by centrifugation. Remove serum and assay immediately or aliquot and store samples at -80°C (preferably) or -20°C. Avoid multiple freeze-thaw cycles
- **Plasma samples** - Blood should be collected into a container with an anticoagulant and then centrifuged. Assay immediately or aliquot and store samples at -80°C (preferably) or -20°C. Avoid multiple freeze-thaw cycles.
- **Urine samples** - Collect mid-stream using a sterile or clean urine collector. Centrifuge to remove cell debris. Assay immediately or aliquot and store samples at -80°C (preferably) or -20°C. Avoid multiple freeze-thaw cycles.
- **Known interfering substances** - Azide and thimerosal at concentrations higher than 0.1% inhibit the enzyme reaction.

ASSAY PROTOCOL

DILUTION OF SAMPLES

The assay for quantification of microalbumin in samples requires that each test sample be diluted before use. All samples should be assayed in **duplicate** each time the assay is performed. The recommended dilutions are only suggestions. Dilutions should be based on the expected concentration of the unknown samples such that the diluted sample falls within the dynamic range of the calibration curve. For absolute quantification of samples that yield results outside the range of the calibration curve, a lesser or greater dilution might be required. **If unsure of sample level, a serial dilution with one or two representative samples before running the entire plate is highly recommended.**

- **Serum and Plasma samples** - Recommended starting dilution is **1:1,000,000**. To prepare a 1:1,000,000 dilution of a sample, transfer 2 µL of sample to 1,998 µL of 1X diluent. This gives you a 1:1,000 dilution. Next dilute the 1:1,000 dilution by transferring 2 µL into 1,998 µL of 1X diluent. This gives you a 1:1,000,000 dilution. Mix thoroughly at each stage.
- **Urine samples**- Recommended starting dilution is **1:500**. To prepare a 1:500 dilution of a sample, transfer 2 µL of sample to 998 µL of 1X diluent. This gives you a 1:500 dilution. Thoroughly mix at each stage

REAGENT PREPARATION

1. **Bring all reagents to room temperature (RT) before use.**
2. **Diluent Concentrate**
The **Diluent** solution supplied is a 5X concentrate and must be diluted **1:5** with distilled or de-ionized water. (1 part buffer concentrate, 4 parts dH₂O)
3. **Wash Solution Concentrate**
The **Wash Solution** supplied is a 20X concentrate and must be diluted **1:20** with distilled or deionized water (1 part buffer concentrate, 19 parts dH₂O). Crystal formation in the concentrate is

not uncommon when storage temperatures are low. Warming of the concentrate to 30-35°C before dilution can dissolve crystals.

4. **Enzyme-Antibody Conjugate Concentrate**
Calculate the required amount of working conjugate solution for each microtiter plate test strip by adding 10 µL **Enzyme-Antibody Conjugate** to 990 µL of **1X Diluent** for each test strip to be used for testing. Dilute immediately before use and protect from light. Mix uniformly, but gently to avoid forming foam.
5. **TMB Substrate Solution**
Ready to use as supplied.
6. **Stop Solution**
Ready to use as supplied.
7. **Microtiter Plate**
Ready to use as supplied. Unseal the microtiter pouch and remove the plate from the pouch. Remove all of the strips and wells that will not be used in the assay and place back in pouch and re-seal along with desiccant.
8. **Rat Microalbumin Calibrator**
Prepare according to the lot specific **Certificate of Analysis**
9. **Positive Control**
The concentration and recommended dilution are provided on the **Control Certificate**. Before use, briefly centrifuge the **Positive Control** to allow all of the liquid to collect in the bottom of the vial.

PROCEDURE

1. **All samples and Calibrators should be assayed in duplicate.**
2. The **Calibrators** and the test sample(s) should be loaded into the ELISA wells as quickly as possible to avoid a shift in optical density (OD) readings. Using a multichannel pipette would reduce this occurrence.

Pipette 100 µL of

Calibrator 0	(0.0 ng/mL) in duplicate
Calibrator 1	(6.25 ng/mL) in duplicate
Calibrator 2	(12.5 ng/mL) in duplicate
Calibrator 3	(25 ng/mL) in duplicate
Calibrator 4	(50 ng/mL) in duplicate
Calibrator 5	(100 ng/mL) in duplicate
Calibrator 6	(200 ng/mL) in duplicate
Calibrator 7	(400 ng/mL) in duplicate

3. Pipette 100 µL of sample (in duplicate) into pre-designated wells.
4. Incubate the **Microtiter Plate** at RT for **thirty** (30 ± 2) minutes. Keep the plate covered and level during incubation.
5. Following incubation, aspirate the contents of the wells.
6. Completely fill each well with appropriately diluted **Wash Solution** and aspirate. Repeat three times, for a total of four washes. If washing manually; completely fill wells with diluted **Wash Solution**, invert the plate and pour/shake out the contents in a waste container. Follow this by sharply striking the wells on absorbent paper to remove residual **Wash Solution**. Repeat 3 times for a total of 4 washes.
7. Pipette 100 µL of appropriately diluted **Enzyme-Antibody Conjugate** to each well. Incubate at RT for **thirty** (30 ± 2) minutes. Keep the plate covered in the dark and level during incubation.
8. Wash and blot the wells as described in Steps 5 and 6.

9. Pipette 100 µL of **TMB Substrate** Solution into each well.
10. Incubate in the **dark** at **RT** for precisely **ten** (10) minutes.
11. After **ten** (10) minutes, add 100 µL of **Stop Solution** to each well.
12. Determine the absorbance (at 450 nm) of the contents of each well within **thirty** (30) minutes.
Calibrate the plate reader to the manufacturer's specifications.

RESULTS

1. Subtract the average background value (average absorbance reading of Calibrator 0) from the test values for each sample.
2. Average the duplicate readings for each calibrator and use the results to construct a Calibration Curve. Construct the curve by reducing the data using computer software capable of generating a four-parameter logistic curve fit. A second order polynomial (quadratic) or other curve fits may also be used; however, they will be a less precise fit of the data.
3. Interpolate test sample values from calibration curve. Correct for sample dilution factor to arrive at the microalbumin concentration in the original sample.

FOR RESEARCH USE ONLY

KAMIYA BIOMEDICAL COMPANY

12779 Gateway Drive, Seattle, WA 98168
Tel: (206) 575-8068 Fax: (206) 575-8094
Email: LifeScience@k-assay.com
www.k-assay.com