

**KAMIYA BIOMEDICAL COMPANY**

# Dog N-Terminal Pro-Brain Natriuretic Peptide (NT- ProBNP) ELISA

**For the quantitative determination of dog NT-ProBNP in  
serum, plasma, tissue homogenates, cell lysates, cell culture  
supernates and other biological fluids.**

**Cat. No. KT-23770**

**For Research Use Only.**

**Product Information**  
**Dog N-Terminal Pro-Brain Natriuretic Peptide (NT-ProBNP) ELISA**  
 Cat. No. KT-23770

## INTENDED USE

The kit is a competitive inhibition enzyme immunoassay technique for the *in vitro* quantitative measurement of dog NT-ProBNP in serum, plasma, tissue homogenates, cell lysates, cell culture supernates and other biological fluids. For research use only.

## COMPONENTS

Reagents	Quantity
Pre-coated, ready to use 96-well strip plate	1
Calibrator	2
Calibrator Diluent	1 × 20 mL
Detection Reagent A	1
Detection Reagent B	1 × 120 µL
Assay Diluent A	1 × 12 mL
Assay Diluent B	1 × 12 mL
TMB Substrate	1 × 9 mL
Stop Solution	1 × 6 mL
Wash Buffer (30X concentrate)	1 × 20 mL
Reagent Diluent	1 × 300 µL
Plate sealer for 96 wells	4

## MATERIALS REQUIRED BUT NOT SUPPLIED

1. Microplate reader with 450 ± 10 nm filter.
2. Single or multi-channel pipettes with high precision and disposable tips.
3. Microcentrifuge Tubes.
4. De-ionized or distilled water.
5. Absorbent paper for blotting the microplate.
6. Container for Wash Solution.
7. 0.01 mol/L (or 1x) Phosphate Buffered Saline (PBS), pH 7.0-7.2.

## STORAGE

All the reagents should be kept according to the labels on vials. The **Calibrator**, **Detection Reagent A**, **Detection Reagent B**, **Reagent Diluent** and the **96-well strip plate** should be stored at -20°C upon receipt while the others should be at 4°C. The unused strips should be kept in a sealed bag with the desiccant provided to minimize exposure to damp air. Opened test kits will remain stable for 1 month, provided it is stored as described above.

## TEST PRINCIPLE

This assay employs the competitive inhibition enzyme immunoassay technique. A monoclonal antibody specific to NT-ProBNP has been pre-coated onto a microplate. A competitive inhibition reaction is launched between biotin labeled NT-ProBNP and unlabeled NT-ProBNP (Calibrators or samples) with the pre-coated antibody specific to NT-ProBNP. After incubation the unbound conjugate is washed off. Next, avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. The amount of bound HRP conjugate is reverse proportional to the concentration of NT-ProBNP in the

sample. After addition of the substrate solution, the intensity of color developed is reverse proportional to the concentration of NT-ProBNP in the sample.

## **SAMPLE COLLECTION AND STORAGE**

### **Serum**

Use a serum separator tube and allow samples to clot for two hours at room temperature or overnight at 4°C before centrifugation for 20 minutes at approximately 1,000 x g. Assay freshly prepared serum immediately or store samples in aliquot at -20°C or -80°C for later use. Avoid repeated freeze/thaw cycles.

### **Plasma**

Collect plasma using EDTA or heparin as an anticoagulant. Centrifuge samples for 15 minutes at 1,000 x g at 4°C within 30 minutes of collection. Remove plasma and assay immediately or store samples in aliquot at -20°C or -80°C for later use. Avoid repeated freeze/thaw cycles.

### **Tissue homogenates**

The preparation of tissue homogenates will vary depending upon tissue type.

1. Tissues were rinsed in ice-cold PBS to remove excess blood thoroughly and weighed before homogenization.
2. Minced the tissues to small pieces and homogenized them in fresh lysis buffer (different lysis buffer needs to be chosen based on subcellular location of the target protein) (w:v = 1:20-1:50, e.g. 1 mL lysis buffer is added in 20-50 mg tissue sample) with a glass homogenizer on ice (Micro Tissue Grinders works, too).
3. The resulting suspension was sonicated with an ultrasonic cell disrupter till the solution is clarified.
4. Then, the homogenates were centrifuged for 5 minutes at 10,000 x g. Collect the supernates and assay immediately or aliquot and store at ≤-20°C.

### **Cell Lysates**

Cells need to be lysed before assaying according to the following directions.

1. Adherent cells should be washed by cold PBS gently, and then detached with trypsin, and collected by centrifugation at 1,000 x g for 5 minutes (suspension cells can be collected by centrifugation directly).
2. Wash cells three times in cold PBS.
3. Resuspend cells in fresh lysis buffer with concentration of 10<sup>7</sup> cells/mL. If it is necessary, the cells could be subjected to ultrasonication till the solution is clarified.
4. Centrifuge at 1,500 x g for 10 minutes at 4°C to remove cellular debris. Assay immediately or aliquot and store at ≤-20°C.

### **Cell culture supernates and other biological fluids**

Centrifuge samples for 20 minutes at 1,000 x g. Collect the supernates and assay immediately or store samples in aliquot at -20°C or -80°C for later use. Avoid repeated freeze/thaw cycles.

### **Note:**


1. Samples to be used within 5 days may be stored at 4°C, otherwise samples must be stored at -20°C (≤1 month) or -80°C (≤2 months) to avoid loss of bioactivity and contamination.
2. When performing the assay, bring samples to room temperature.
3. Sample hemolysis will influence the result, so hemolytic specimen should not be used.
4. It is highly recommended to use serum instead of plasma for the detection based on quantity of our in-house data.

## **REAGENT PREPARATION**

Bring all kit components and samples to room temperature (18-25°C) before use. If the kit will not be used up in one time, please only take out strips and reagents for present experiment, and leave the remaining strips and reagents in required condition.

**Calibrator**

Reconstitute the **Calibrator** with 1.0 mL of **Calibrator Diluent**, kept for 10 minutes at room temperature, shake gently (not to foam). The concentration of the calibrator in the stock solution is 50,000 pg/mL. Please firstly dilute the stock solution to 5,000 pg/mL and the diluted calibrator serves as the highest calibrator (5,000 pg/mL). Then prepare 5 tubes containing 0.6 mL Calibrator Diluent and produce a triple dilution series according to the picture shown below. Mix each tube thoroughly before the next transfer. Set up 5 points of diluted calibrator such as 5,000 pg/mL, 1,666.7 pg/mL, 555.6 pg/mL, 185.2 pg/mL, 61.7 pg/mL, and the last EP tubes with **Calibrator Diluent** is the blank as 0 pg/mL.

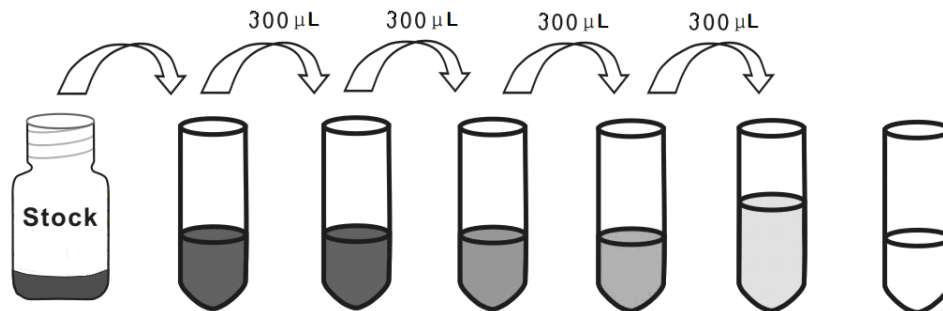

<b>Tube</b>	<b>1</b>	<b>2</b>	<b>3</b>	<b>4</b>	<b>5</b>	<b>6</b>	<b>7</b>
<b>pg/mL</b>	<b>50,000</b>	<b>5,000</b>	<b>1,666.7</b>	<b>555.6</b>	<b>185.2</b>	<b>61.7</b>	<b>0</b>

**Detection Reagent A**

Reconstitute the **Detection Reagent A** with 150 µL of **Reagent Diluent**, kept for 10 minutes at room temperature, shake gently (not to foam). Dilute to the working concentration with **Assay Diluent A** (1:100).

**Detection Reagent B**

Briefly spin or centrifuge the stock Detection B before use. Dilute to the working concentration with **Assay Diluent B** (1:100).

**Wash Solution**

Dilute 20 mL of Wash Solution Concentrate (30X) with 580 mL of de-ionized or distilled water to prepare 600 mL of Wash Solution (1X).

**TMB substrate**

Aspirate the needed dosage of the solution with sterilized tips and do not dump the residual solution into the vial again.

**Note:**

1. Making serial dilution in the wells directly is not permitted.
2. Prepare calibrator within 15 minutes before assay. Please do not dissolve the reagents at 37°C directly.
3. Please carefully reconstitute Calibrators or working Detection Reagent A and B according to the instruction, and avoid foaming and mix gently until the crystals are completely dissolved. To minimize imprecision caused by pipetting, use small volumes and ensure that pipettors are calibrated. It is recommended to suck more than 10 µL for one pipetting.
4. The reconstituted Calibrators, Detection Reagent A and Detection Reagent B can be **used only once**.
5. If crystals have formed in the Wash Solution concentrate (30X), warm to room temperature and mix gently until the crystals are completely dissolved.
6. Contaminated water or container for reagent preparation will influence the detection result.
7. Detection Reagent A and B are sticky solutions, therefore, slowly pipette them to reduce the volume errors.

## SAMPLE PREPARATION

1. We are only responsible for the kit itself, but not for the samples consumed during the assay. The user should calculate the possible amount of the samples used in the whole test. Please reserve sufficient samples in advance.
2. Please predict the concentration before assaying. If values for these are not within the range of the calibration curve, users must determine the optimal sample dilutions for their particular experiments. Sample should be diluted by PBS.
3. If the samples are not indicated in the manual, a preliminary experiment to determine the validity of the kit is necessary.
4. Tissue or cell extraction samples prepared by chemical lysis buffer may cause unexpected ELISA results due to the impacts from certain chemicals.
5. Due to the possibility of mismatching between antigen from other origin and antibody used in our kits (e.g., antibody targets conformational epitope rather than linear epitope), some native or recombinant proteins from other manufacturers may not be recognized by our products.
6. Influenced by the factors including cell viability, cell number or sampling time, samples from cell culture supernates may not be detected by the kit.
7. Fresh samples without long time storage are recommended for the test. Otherwise, protein degradation and denaturalization may occur in those samples and finally lead to wrong results.

## ASSAY PROCEDURE

1. Determine wells for diluted calibrator, blank and sample. Prepare 5 wells for calibrator points, 1 well for blank. Add 50  $\mu$ L each of dilutions of calibrator (read Reagent Preparation), blank and samples into the appropriate wells, respectively. And then add 50  $\mu$ L of Detection Reagent A to each well immediately. Shake the plate gently (using a microplate shaker is recommended). Cover with a Plate sealer. Incubate for 1 hour at 37°C. Detection Reagent A may appear cloudy. Warm to room temperature and mix gently until solution appears uniform.
2. Aspirate the solution and wash with 350  $\mu$ L of 1X Wash Solution to each well using a squirt bottle, multi-channel pipette, manifold dispenser or autowasher, and let it sit for 1~2 minutes. Remove the remaining liquid from all wells completely by snapping the plate onto absorbent paper. Repeat 3 times. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against absorbent paper.
3. Add 100  $\mu$ L of Detection Reagent B working solution to each well. Incubate for 30 minutes at 37°C after covering it with the Plate sealer.
4. Repeat the aspiration/wash process for total 5 times as conducted in step 2.
5. Add 90  $\mu$ L of Substrate Solution to each well. Cover with a new Plate sealer. Incubate for 10 - 20 minutes at 37°C (Don't exceed 30 minutes). Protect from light. The liquid will turn blue by the addition of Substrate Solution.
6. Add 50  $\mu$ L of Stop Solution to each well. The liquid will turn yellow by the addition of Stop solution. Mix the liquid by tapping the side of the plate. If color change does not appear uniform, gently tap the plate to ensure thorough mixing.
7. Remove any drop of water and fingerprint on the bottom of the plate and confirm there is no bubble on the surface of the liquid. Then, run the microplate reader and conduct measurement at 450 nm immediately.

**Note:**

1. **Assay preparation:** Keep appropriate numbers of wells for each experiment and remove extra wells from microplate. Rest wells should be resealed and stored at -20°C.
2. **Samples or reagents addition: Please use the freshly prepared Calibrator.** Please carefully add samples to wells and mix gently to avoid foaming. Do not touch the well wall. For each step in the procedure, total dispensing time for addition of reagents or samples to the assay plate **should not exceed 10 minutes**. This will ensure equal elapsed time for each pipetting step, without interruption. Duplication of all calibrators and specimens, although not required, is recommended. To avoid cross-contamination, change pipette tips between additions of calibrators, samples, and reagents. Also, use separated reservoirs for each reagent.
3. **Incubation:** To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary. Do not allow wells to sit uncovered for extended periods between incubation steps. Once reagents are added to the well strips, DO NOT let the strips DRY at any time during the assay. Incubation time and temperature must be controlled.
4. **Washing:** The wash procedure is critical. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Solution by aspirating or decanting and remove any drop of water and fingerprint on the bottom of the plate. Insufficient washing will result in poor precision and falsely elevated absorbance reading.
5. **Controlling of reaction time:** Observe the change of color after adding **TMB Substrate** (e.g. observation once every 10 minutes), if the color is too deep, add **Stop Solution** in advance to avoid excessively strong reaction which will result in inaccurate absorbance reading.
6. **TMB Substrate** is easily contaminated. Please protect it from light.
7. The environment humidity which is less than 60% might have some effects on the final performance, therefore, a humidifier is recommended to be used at that condition.

## CALCULATION OF RESULTS

This assay employs the competitive inhibition enzyme immunoassay technique, so there is an inverse correlation between NT-ProBNP concentration in the sample and the assay signal intensity.

Average the duplicate readings for each calibrator, control, and samples. Create a calibration curve with the log of NT-ProBNP concentration on the y-axis and absorbance on the x-axis. Draw a best fit curve through the points and it can be determined by regression analysis. If samples have been diluted, the concentration read from the calibration curve must be multiplied by the dilution factor.

## PERFORMANCE

### Detection Range

61.7- 5,000 pg/mL.

The calibration curve concentrations used for the ELISA's were 5,000 pg/mL, 1,666.7 pg/mL, 555.6 pg/mL, 185.2 pg/mL, 61.7 pg/mL.

### Sensitivity

The minimum detectable dose of NT-ProBNP is typically less than 23.5 pg/mL.

The sensitivity of this assay, or Lower Limit of Detection (LLD) was defined as the lowest protein concentration that could be differentiated from zero. It was determined by subtracting two standard deviations to the mean optical density value of twenty zero calibrator replicates and calculating the corresponding concentration.

### Specificity

This assay has high sensitivity and excellent specificity for detection of NT-ProBNP.

No significant cross-reactivity or interference between NT-ProBNP and analogues was observed.

**Note:**

Limited by current skills and knowledge, it is impossible for us to complete the cross-reactivity detection between NT-ProBNP and all the analogues, therefore, cross reaction may still exist.

**ASSAY PROCEDURE SUMMARY**

1. Prepare all reagents, samples and calibrators;
2. Add 50 µL calibrator or sample to each well. And then add 50 µL prepared Detection Reagent A immediately. Shake and mix. Incubate 1 hour at 37°C;
3. Aspirate and wash 3 times;
4. Add 100 µL prepared Detection Reagent B. Incubate 30 minutes at 37°C;
5. Aspirate and wash 5 times;
6. Add 90 µL Substrate Solution. Incubate 10-20 minutes at 37°C;
7. Add 50 µL Stop Solution. Read at 450 nm immediately.

**IMPORTANT NOTE**

1. The final experimental results will be closely related to validity of the products, so the kit should be used prior to the expiration date. And please store the kits exactly according to the instruction.
2. Kits from different batches may be a little different in detection range, sensitivity and color developing time. Please perform the experiment exactly according to the instruction attached in kit while electronic ones from our website is only for reference.
3. Do not mix or substitute reagents from one kit lot to another. Use only the reagents supplied by manufacturer.
4. Protect all reagents from strong light during storage and incubation. All the bottle caps of reagents should be covered tightly to prevent the evaporation and contamination of microorganism. TMB Substrate should remain colorless till it is reacted with the enzyme which binds to the microplate.
5. There may be some foggy substance in the wells when the plate is opened at the first time. It will not have any effect on the final assay results. Do not remove microplate from the storage bag until needed.
6. Wrong operations during the reagents preparation and loading, as well as incorrect parameter setting for the plate reader may lead to incorrect results. A microplate reader with a bandwidth of 10 nm or less and an optical density range of 0-3 O.D. at  $450 \pm 10$  nm wavelength is acceptable for use in absorbance measurement. Please read the instruction carefully and adjust the instrument prior to the experiment.
7. Variation in sample preparation and each step of experimental operation may cause different results. In order to get better reproducible results, the operation of each step in the assay should be controlled.
8. Each kit has been strictly passed Q.C test. However, results from end users might be inconsistent with our in-house data due to some unexpected transportation conditions or different lab equipments. Intra-assay variance among kits from different batches might arise from above factors, too.
9. Kits from different manufacturers with the same item might produce different results, since we haven't compared our products with other manufacturers.
10. The calibrator of the kit and immunogen used for antibody preparation are commonly recombinant proteins, as different fragments, expression systems, purification methods might be used in recombinant protein preparation, we can not guarantee the kit could detect recombinant protein from other companies. So, it is not recommended to use the kit for the detection of recombinant protein.
11. Please predict the concentration of target molecules in samples, or arrange a preliminary experiment, it is a good way to solve specific problem, e.g. the concentration of samples are beyond the detection range of the kit.
12. The kit might not be suitable for detection of samples from some special experiment, for instance, knock-out experiments, due to their uncertainty of effectiveness.
13. The Stop Solution suggested for use with this kit is an acid solution. Wear eye, hand, face, and clothing protection when using this material.

**FOR RESEARCH USE ONLY.**

**KAMIYA BIOMEDICAL COMPANY**

12779 Gateway Drive, Tukwila, WA 98168  
Tel: (206) 575-8068 Fax: (206) 575-8094  
Email: LifeScience@k-assay.com  
[www.k-assay.com](http://www.k-assay.com)